

HOW TO OBTAIN SUPPLIES

CATHETER DESCRIPTION	FP10 CODE	PIP CODE
HydroSil male 10Ch	63610	301-4396
HydroSil male 12Ch	63612	301-4404
HydroSil male 14Ch	63614	301-4412
HydroSil male 16Ch	63616	301-4420
HydroSil male 18Ch	63618	301-4438

For fast, reliable delivery of HydroSil male catheters, please call the Script-easy confidential prescription service.

Script-easy 
CONFIDENTIAL PRESCRIPTION DISPENSING SERVICE

STEP 1 Patient Registration

Patient registration by phone, fax or post



STEP 2 Script-easy Service

The Script-easy service liaises with patient and GP to obtain prescription



STEP 3 Home Delivery

Home delivery by confidential carrier every month



Contact Script-easy by Freephone on 0800 0121699



DISTRIBUTED BY
Rochester Medical Ltd
10 Commerce Way
Lancing, West Sussex BN15 8TA

Freephone: 0800 0322755
Phone: +44 (0)1903 875055
Fax: +44 (0)1903 875085

MANUFACTURED BY
Rochester Medical Corporation, Stewartville, MN 55976 USA

© 2007 Rochester Medical Corporation. All rights reserved. HydroSil and Script-easy are trademarks of Rochester Medical Corporation. Rochester Medical is a registered trademark of Rochester Medical Corporation. Patents pending.

540197

HydroSil[®]
male

Male Hydrophilic Intermittent Catheter

More Comfortable.
More Convenient.
More Confident.



HydroSil[®] male

Male Hydrophilic Intermittent Catheter

HydroSil has a unique hydrophilic coating that provides a soft cushion of water for virtually friction-free insertion and removal and offers the ideal qualities for optimal feel and comfort.

Four Comfort Sized Eyes

Enables optimal drainage

Smooth, Tapered, Seamless Tip

Atraumatic insertion

Immediate Hydration

No uncertainty

Silicone Construction

Toxin-free and environmentally friendly


Pure, Sterile Water

Eliminates concerns about clean water source

The Next Generation in Catheter Comfort and Design

Core Technology

Three layer construction designed to provide optimal handling characteristics while remaining comfortable to use

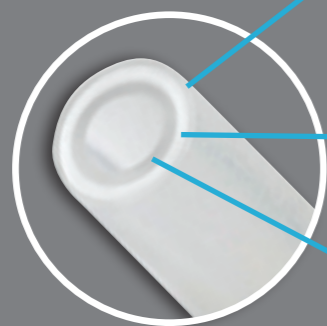

1: Outer Surface Area

Soft, smooth silicone which is gentle on sensitive urethral tissue providing atraumatic insertion and withdrawal

2: Center Core

Firm for optimal handling characteristics and urethral navigation

3: Inner Layer

Hydrophobic surface designed to maximise drainage flow and allow the catheter to more freely navigate the urethra

Nurse HydroSil “take-home” packs include 15 catheters, the HydroSil Users’ Guide, and details of the Script-*easy* prescription service

Patient HydroSil prescriptions through Script-*easy* come with thoughtful accessories, including a catheter carrying pouch, wet wipes and a non-alcohol based hand sanitiser

